


- 3 section treatment platform
- Electric hi/lo from 46 to 94 cm controlled by foot pedals accessible from all sides of the table
- Gas-assisted head section pivots up and down independent of the chest section
- Gas-assist elevates the foot section with the quick-release grab bar. The leg cushion can be raised to any angle desired up to 80 degrees
- Patient positioning system with i-skin integral Skin foam standard
- Dual-lever retractable caster system
- Engage the optional, motorized postural drain feature with the push of a button
- Easily convert to Trendelenberg position with a quick raise of the gas-assisted grab bar



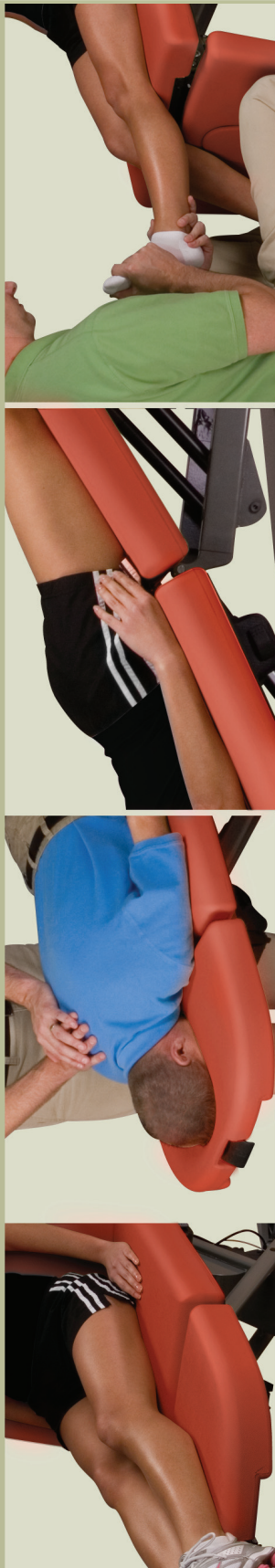

ORDERING INFORMATION

Region 3

3330 – Region 3 treatment platform

3331 – Region 3D treatment platform
With postural drainage

TECHNICAL SPECIFICATIONS

Mains Power: 100-240 V, 50-60 Hz

Standby Consumption: 7 W

Input: 230 V / 1 A / 50 Hz, 120 V / 1.7 A / 60 Hz

Weight: 141 kg (310 lbs.)

Dimensions (LxWxH): 200x71 x 45.7–94 cm
(79" x 28" x 18"-37")

Frame Shipping Dimensions (LxWxH):
152 x 74 x 62.5cm (60" x 29" x 25")

Cushion Shipping Dimensions (LxWxH):
69 x 16 x 108cm (27" x 6" x 42.5")

Lifting Capacity: 204 kg (450 lbs)

Electrical Safety Class: Class II, Type B

Safety Tests: IEC/UL/EN 60601-1, IEC
60601-1-1, IEC 60601-1-2,
CAN/USA C22.2 No. 601.1-M90 w/ A2

SKIN COLORS

Sand

Sage

Sunflower

Terra Cotta

Pacific Blue

Maroon

Chocolate

Midnight

TABLE ARTICULATIONS

TABLE SCHEMATIC

6